




Operating Instructions

863 Valley View Road, Eighty Four, PA 15330
Tel: 724-941-9701 • www.skcinc.com

OVS Tube Holder Cat. No. 224-29V

Description

The SKC OVS Tube Holder prevents accidental breakage of the sample tube during use. It is designed for use with the OSHA Versatile Sampler (OVS) tubes (SKC Cat. Nos. 226-30-16, 226-56, 226-57, 226-58, 226-99). The OVS Tube Holder consists of a standard tube holder fitted with a special durable protective cover for the OVS tubes and a 3-foot (0.9-meter) length of Tygon tubing. A collar clip fastened to the protective cover allows for easy attachment to the worker. A flexible metal coil surrounds the tubing just below the holder, and allows the tubing to bend without blocking air flow.

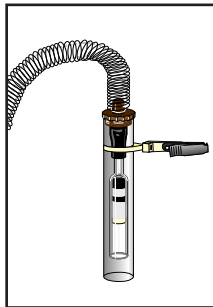

How to Use

To use the OVS Tube Holder, unscrew the cylindrical white protective cover from the dark brown tube holder. Insert the small end of the OVS tube securely into the rubber sleeve of the tube holder. Replace the protective cover. Attach the free end of Tygon tubing to the inlet of the sampling pump. Clip the protective cover to the clothing in the worker's breathing zone with the inlet facing down. Sample at the recommended flow rate.